



Dynamic screen with image and sound

GEWA One – remote control with telephone function

GEWA One is a remote control that helps people to manage equipment in their own home, such as the TV, music system, lamps, profiling beds, door and window openers.

GEWA One also works as a mobile phone enabling the user to keep in touch with friends, family and be able to quickly reach a carer through either a call or text message.

The user can choose whether to control the device through either the touch screen or navigate with an external switch. Symbols, colours and voice recordings help the person quickly find the right button.

Art.nr 4616301



ABILIA

A COMPLETE SOLUTION


GEWA One is a remote control with telephone function, which gives the user the opportunity to call and text. You can customize how pages and buttons should look and have the opportunity to add speed dials to contacts for quicker access.

CONTROL IT YOUR WAY

You only need to tap the screen lightly to activate a button. You can utilise larger buttons for a larger target area, instead of having to navigate the remote with an external control switch.

EASY TO SEE, HEAR AND UNDERSTAND

Symbols, different colours and recorded speech on the buttons can be programmed for easier navigation and choose the number of buttons that a page should display.

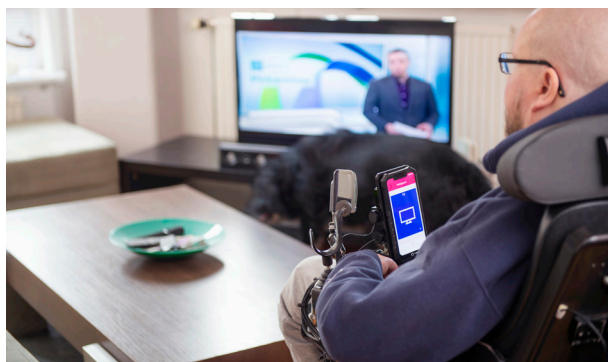

SMART TECHNOLOGY

GEWA One can record signals from a TV's IR remote control. GEWA One can also transmit signals to a GEWA bed receiver for bed control and coded signals for a front door or external window. Insert a SIM card for calling and texting with loved ones.

SUPPORT AND SECURITY

When a technician is going to customize GEWA One for a user, this is done smoothly in our myAbilia web portal. Within myAbilia there is a backup copy of the users remote control setup making some support functions accessed remotely without the need for an engineer to visit.

EASY TO HANDLE

Accessories such as a holder specially designed for the GEWA One allow you to see and manage the remote control from both wheelchair and bed.

COMES WITH

- Mobile unit, restricted to GEWA One
- GEWA Multibox
- GEWA smart cable (Mobile unit-Multibox)
- Charger GEWA One
- USB-cable USB-C to USB-A
- User manual
- Voucher with license key

TECHNICAL INFORMATION GEWA ONE 4616301

MOBILE UNIT		MULTIBOX		ACCESSORIES	ART NR
Screen size	6.6"	Size	108 x 67 x 22 mm	Holder	461198
Size	174.45 x 81.8 x 13.6 mm	Weight	153 g	Mini-arm	451030
Weight	310 g	Connections:	3.5 mm, Micro-USB	Universal arm	621966
OS	Android 12	Control signals:	IR, Radio		

GEWA One is a medical device intended for indoor use.

